

# Bimodal Hearing Aid Fitting Guidelines

WHITEPAPER  
2016

## ABSTRACT

In a bimodal fitting, one ear is stimulated acoustically with a hearing aid and the other is stimulated electrically with a cochlear implant.

To bring bimodal benefits to all children and adults with unilateral implants and audible hearing in the contralateral ear, Oticon has implemented a bimodal fitting guide in the Genie fitting software. As developed by Carisa Reyes, Staff Audiologist at Boys Town National Research Hospital, the bimodal fitting flowchart serves as a guide to clinical audiologists as they navigate the bimodal fitting process. The goal is to provide a logic- and evidence-based method for decision-making, yet keeping in mind the constraints in everyday clinical practice.

Based on the latest knowledge on bimodal research, this paper explains the rationale, the recommended strategies, the procedures and the caveats of the bimodal fitting.

- Candidates for bimodal fitting
- When to fit the hearing aid
- How to fit the hearing aid for bimodal patients
- Bimodal flowchart
- Evaluation of benefits
- Case Studies

*Acknowledgments: Thank you to Marc A. Brennan, Ph.D. and Lisa S. Davidson, Ph.D. for reviewing the manuscript.*



*The bimodal flowchart takes into account the four fitting approaches: wideband fitting, restricted bandwidth fitting, use of frequency lowering and loudness balancing. The chart is available in Genie when you fit Oticon Dynamo or Sensei Super Power.*



**BOYS TOWN**  
National Research  
Hospital



Carisa Reyes, Au.D., CCC-A.  
Staff Audiologist at Boys Town  
National Research Hospital

**oticon**  
PEOPLE FIRST

### Candidates for bimodal fitting:

More than half of those receiving a cochlear implant have aidable hearing in the non-implanted ear<sup>1</sup>. If these recipients are fitted with a hearing aid in the non-implanted ear, access to bilateral and binaural cues (such as those arising from head shadow and redundancy) as well as complementary acoustic cues (such as fundamental frequency) may enhance cochlear implant performance. In this guideline, bimodal stimulation refers to the use of a cochlear implant in one ear and a hearing aid in the opposite ear and bimodal benefit refers to performance improvement with the hearing aid over performance with the cochlear implant alone. Even recipients with significant hearing loss in the non-implanted ear demonstrate bimodal benefit<sup>2-4</sup>. For those who do not obtain bimodal benefit for speech recognition in noise<sup>5</sup>, other benefits such as enhanced music and pitch perception may still occur<sup>6-9</sup>. In addition, many studies have documented more natural sound quality and improved ease of listening with bimodal stimulation<sup>3,10-15</sup>. Finally, bimodal stimulation can facilitate spoken language and literary skills in young children who subsequently receive a second side cochlear implant<sup>16-18</sup>. The provision of bilateral hearing is considered the standard of care for cochlear implant recipients<sup>19,20</sup>. Therefore, all unilateral cochlear implant recipients who have some degree of aidable residual hearing in the contralateral ear should be considered candidates for bimodal stimulation.

### *Fitting a hearing aid in the non-implanted ear gives acoustic cues that may enhance CI performance*

### When to fit the hearing aid:

The recommended time between activation of the cochlear implant and bimodal stimulation varies for different clinics<sup>21,22</sup>. Some argue that delaying bimodal stimulation can facilitate the adjustment to the cochlear implant. However, cessation of hearing aid use or delaying the hearing aid fitting has not been shown to improve *long-term* outcomes. Many cochlear implant candidates are already fitted with bilateral amplification so continued use of the contralateral hearing aid after implantation will allow for continued bilateral stimulation. Immediate hearing aid use will allow the user to obtain the benefits associated with bimodal stimulation sooner rather than later. Continued use of amplification may also facilitate the patient's adjustment to the cochlear implant. There may be isolated cases where this approach would not be appropriate, but in general, bimodal fitting is recommended as soon as possible.

If further optimization of the bimodal fitting is to be undertaken, it has been recommended that this be completed once the cochlear implant program is stable<sup>23,24</sup>. However, clinical judgment may dictate that this be completed sooner (as long as the cochlear implant is set at a comfortable level) because it may take several months to arrive at a stable program.

### How to fit the hearing aid for bimodal patients:

Unfortunately, studies have found that many bimodal users have hearing aids that are fit sub-optimally (i.e., set below targets)<sup>25,26</sup>. It is beyond the scope of this paper to outline the standard hearing aid fitting process, but at a minimum, recommended guidelines for proper selection and verification of amplification should be followed to ensure maximum audibility and comfort of sounds of varying input levels<sup>27,28</sup>.


### *Not a "one-size-fits-all" procedure. The approach that provides the greatest benefit will likely vary among patients*

It has been suggested that optimizing the hearing aid fitting for bimodal use (by considering frequency response and loudness balance) may result in greater benefit over standard hearing aid fitting<sup>2,3,29,30</sup>. However, there is currently no single approach to bimodal fitting that is universally accepted<sup>21,22</sup>. Additionally, there is conflicting evidence that indicates that further optimization of the hearing aid for bimodal use beyond fitting to target using proper verification methods does not result in superior outcomes<sup>31</sup>. Nevertheless, certain fitting strategies may be considered and explored for bimodal recipients as the literature suggests these approaches have the potential for improving outcomes. It should be noted that the approach that provides the greatest benefit or that subjectively provides the most satisfactory sound will likely vary among recipients<sup>32</sup>.


### Bimodal decision-making flowchart:

The accompanying Bimodal Fitting Flowchart was developed to serve as a guide for the bimodal hearing aid fitting process. This is not intended to be a "one-size-fits-all" procedure that is applicable to every bimodal listener. The cochlear implant population is extremely diverse and clinical judgment may warrant that a different approach be taken for a particular patient or patient population. For

### Bimodal flowchart.

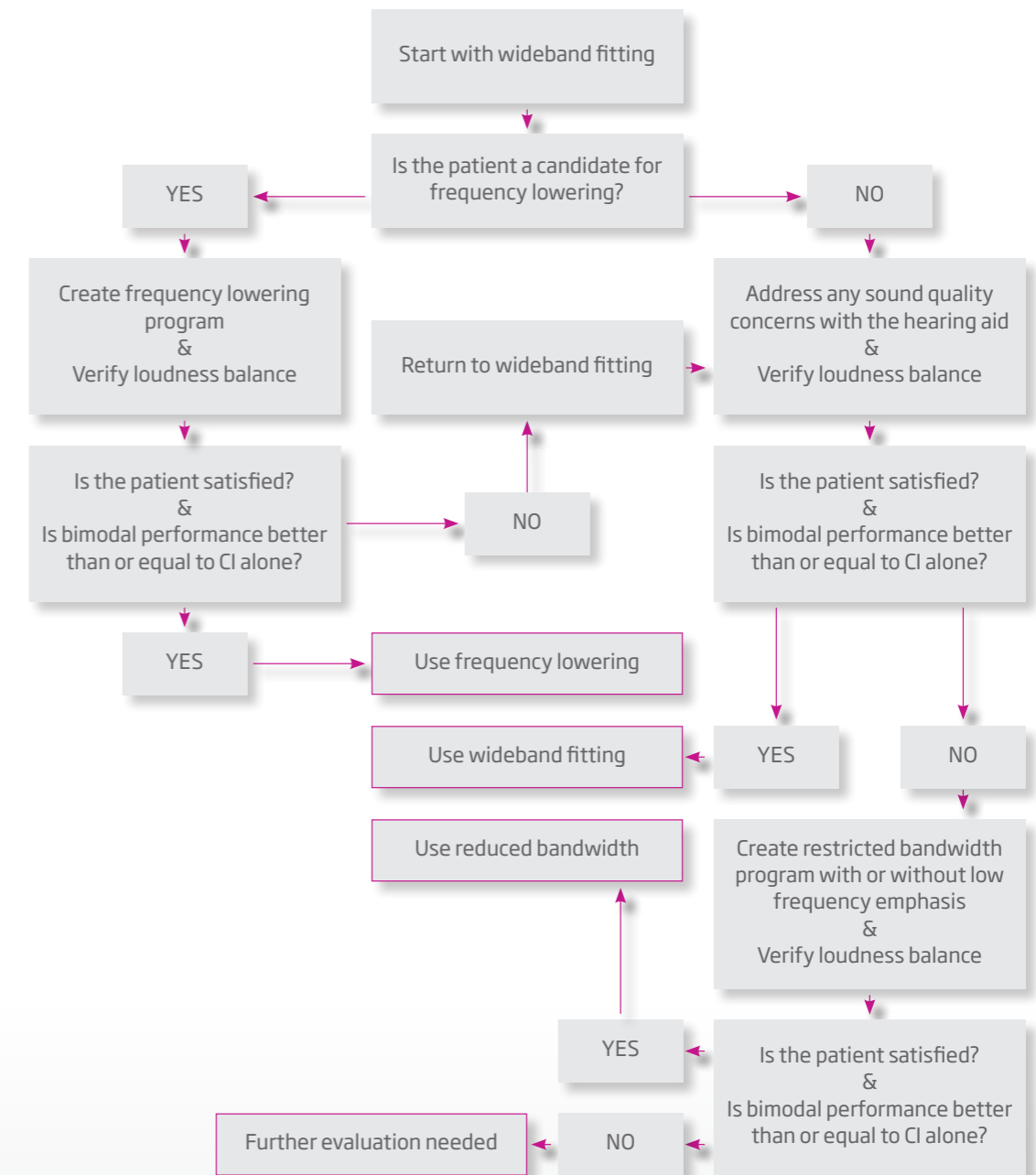

All unilateral cochlear implant recipients with aidable residual hearing in the other ear are candidates for hearing aid use. This flowchart provides an evidence-based, yet practical, method for fitting a hearing aid on a bimodal patient. The flowchart takes into account wideband fitting, restricted bandwidth fitting, use of frequency lowering and loudness balancing.

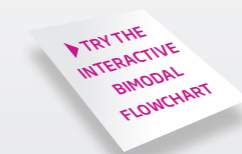

example, there may be varying considerations with adult versus paediatric hearing aid fittings, where in high frequency audibility may be more critical for the paediatric population<sup>33</sup>.

### *A logical and evidence-based method for decision-making, keeping in mind constraints in everyday clinical practice*

The goal of this model is to discuss the various approaches to hearing aid fitting for bimodal users that can be followed using a logical and evidence-based method for decision-making, yet keeping in mind the constraints that can be found in everyday clinical practice (e.g. time constraints, different professionals managing the hearing aid and cochlear implant, etc.). This guideline focuses on the hearing aid fitting. There will be instances where adjustment of the cochlear implant would be more appropriate.

The flowchart incorporates both the hearing aid frequency response and loudness balancing. These are discussed in greater detail below.

#### **Wideband Fitting:**

It is recommended to start with wideband fitting as bimodal benefits have been consistently demonstrated using this standard approach. Prior to obtaining the cochlear implant, a majority of candidates will have amplification already fitted using this approach.

**Rationale:** Apart from optimizing audibility and ensuring listening comfort, the goal of hearing aid fitting for cochlear implant recipients is to allow access to as many potential bilateral, binaural and complementary acoustic cues as possible in order to maximize bimodal benefits. These cues may include high frequency information that can potentially provide inter-aural level difference cues<sup>34,35</sup>. In addition, low frequency acoustic information may provide voice and musical pitch cues that are not transmitted well by the cochlear implant<sup>7,36-38</sup>. Altogether, these additional cues may contribute to improved localization, music perception and speech recognition (particularly in noise). Several studies have demonstrated significant bimodal benefits with this approach versus restricting amplification to the lower frequencies<sup>32,39,40</sup>.

**Recommended Strategy:** Match targets for as wide a bandwidth as possible according to the appropriate prescriptive formula using real ear or simulated real ear measurements.

#### **Restricted Bandwidth with or without low frequency emphasis fitting:**

While it may be prudent to start with wideband fitting, it should be noted that many cochlear implant recipients have significant high frequency hearing losses and/or suspected dead regions in the non-implanted ear. It may not be possible to provide high frequency amplification or to utilize frequency lowering. Even if amplification of the mid to higher frequencies is possible, this may degrade performance in certain patients. Potential advantages over wideband fitting include improved battery life, preventing "off-frequency" listening and feedback reduction.

**Rationale:** Several studies have suggested that the majority of the bimodal benefit for speech perception is obtained from the lower frequencies<sup>4,40-43</sup>. One adult study found that their subjects demonstrated greater bimodal benefit when amplification was not provided beyond the edge frequency of the dead region<sup>44</sup>. More research is needed in this area but the approach of limiting high frequency amplification could be investigated in cases of "bimodal decrement" (poorer performance in the bimodal condition vs. with the cochlear implant alone) or where there is lack of objective and subjective bimodal benefit. This is supported by work by Mok et al. (2006; 2010)<sup>45,46</sup> who found that greater bimodal benefit was found in subjects with poorer mid- and/or high-frequency aided thresholds. Finally, a few authors have suggested or demonstrated that some bimodal users prefer and/or derive benefit from alternative frequency responses, including one that provides additional low frequency emphasis while de-emphasizing the higher frequencies<sup>2,3,24,29</sup>.

**Caveat:** Sound localization may be poorer when utilizing this approach versus with wideband fitting or frequency lowering<sup>32,35</sup>. Because cutting out high frequency amplification will affect access to inter-aural level difference cues, this approach is not recommended as the initial fitting strategy for bimodal users.

**Recommended Strategies:** It should be noted that there is no single accepted method for determining how to restrict high frequency amplification for bimodal users. The following can be considered:


Based on the work of Zhang et al. (2014), the edge frequency of a cochlear dead region is determined using the Sweeping Psychophysical Tuning Curve (SWPTC) Test or the Threshold-Equalizing Noise (TEN) Test. Amplification is only provided up to the edge frequency of the dead region.

The SWPTC or TEN Test may not be clinically available and/or some patients (e.g. young children) may not be able to complete this type of testing. Because studies have shown that dead regions are oftentimes present when thresholds are in the severe to profound range<sup>47</sup>, the dead region could potentially be estimated based on pure tone thresholds. That is, amplification is only provided up to frequencies where thresholds are equal to or better than 80-90 dB HL. Davidson et al. (2015) looked at the older paediatric to young adult population and the cut-off frequency was defined as the lowest frequency where the threshold was > 90 dB HL and the root mean square (RMS) average of the aided speech map for an input of 60 dB SPL fell below that threshold. Zhang et al. (2014), on the other hand, looked at adult subjects and selected the cut-off frequency based on 80 dB HL thresholds and compared this to TEN and SWPTC Test results. They found that cut-offs based on 80 dB thresholds matched with TEN Test results in 5 out of 11 cases and were always higher than cut-offs obtained from SWPTC Test results. This suggests that one may need to experiment with different cut-off frequencies to determine the one that results in the best possible outcomes.

#### **Frequency Lowering:**

In order to try to match frequency bandwidth for contralateral acoustic and electric hearing, high frequency audibility on the hearing aid needs to be maximized. Frequency lowering is now available in many commercial devices and has been used as a strategy to improve high frequency audibility with concurrent improvements in outcomes.

**Rationale:** For hearing aid users, frequency lowering improves high frequency detection and speech recognition<sup>48,49</sup>. For bimodal users, most research has not demonstrated significant benefit over standard fitting for subjects using frequency transposition<sup>50</sup> or frequency compression devices<sup>51-53</sup>. However, most of these studies indicated a high acceptance rate for frequency lowering and no decrement in performance. A recent study looking at older children and young adults<sup>32</sup> points to success using frequency lowering, with many subjects performing best in terms of localization, talker recognition and speech recognition with this approach. In addition, most subjects preferred frequency lowering over wideband and restricted bandwidth fittings.

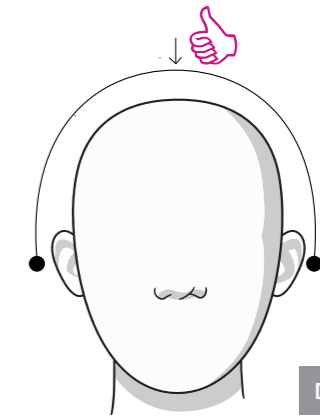

[Download PDF](#)

*Loudness balancing is an attempt to balance loudness between the CI and the hearing aid such as they are judged to be equally loud. The balance point for the electric and acoustic sound is found by increasing and decreasing the level of the hearing aid. Print and attach this picture to a wall. Ask the patient to indicate on the arc where sound is coming from. The balance is achieved when the overall sound produces a sensation that the stimulus is perceived directly in front of the head (arrow). For some patients this procedure may be easier than judgments of the relative loudness of each device.<sup>59</sup>*

**Candidacy:** The Clinical Practice Guidelines for Paediatric Amplification<sup>28</sup> (American Academy of Audiology, 2013) indicate that frequency lowering can be considered if high frequency audibility is not possible using conventional amplification. The potential for improvement of audibility will depend on the patient's hearing loss, type of frequency lowering available and settings chosen. In general, a patient is a candidate for frequency lowering if the hearing aid fitting shows enough audible bandwidth to be able to take the inaudible high frequency components to a frequency where they can be made audible without causing harmful distortion.

**Recommended Strategy:** McCreery et al. (2013) showed that frequency lowering fittings that maximized audible bandwidth resulted in better outcomes versus fittings that did not<sup>49</sup>. The goal of maximizing audibility is also applicable to all forms of frequency lowering.

**Procedure:** Determine Maximum Audible Output Frequency (i.e., the highest frequency in the amplified speech signal that is audible) via real ear or simulated real ear measurements. Then, find the appropriate frequency lowering settings that maximize audibility and minimize distortion by using available tools (e.g. Verifit or the Frequency Lowering Fitting Assistants<sup>48,54</sup>).



### Loudness Balancing:

With properly fitted hearing aids and cochlear implants, sound should be as audible as possible and comfortable for varying input levels in both ears. Occasionally, however, recipients may require adjustments to address loudness concerns. For example, a low frequency gain increase may be needed to compensate for the reduction in perceived loudness when a restricted bandwidth fitting is utilized<sup>32</sup>. Alternately, a gain decrease may be needed due to binaural loudness summation<sup>55</sup>.

Given differences in loudness growth with acoustic versus electric hearing<sup>55</sup> as well as other factors related to asymmetries created by signal processing differences between devices (amplitude compression characteristics, frequency range, noise reduction, etc.)<sup>56,57</sup>, it may be difficult to match loudness for acoustic and electric hearing. Nevertheless, several authors have suggested or demonstrated that greater bimodal benefits for sound localization and speech recognition could be obtained with loudness-balanced devices<sup>2,3,23,55</sup>.

**Recommended Strategy:** There are many procedures for loudness balancing described in the literature<sup>2,23,45,58</sup>. However, these typically require calibrated stimuli, tend to be time-intensive and may be difficult for patients to perform. Dorman et al. (2014)<sup>59</sup> found that although bimodal benefits can be obtained with unequal loudness between ears, the greatest benefit was seen when the acoustic signal was judged as equally loud to just noticeably softer than the CI signal or if the acoustic signal was perceived at “most comfortable loudness” level. The authors concluded that a high degree of precision is not necessary when attempting to balance loudness between devices. Although the procedure described below focuses on adjusting the hearing aid, there may be instances when it is more appropriate to adjust the overall loudness of the cochlear implant.

**Procedure:** If the hearing aid needs to be adjusted, apply gain and/or compression adjustments until a perception of equal (or close to equal) loudness with the cochlear implant is obtained. Alternately, adjustments can also be made when listening with the hearing aid alone so that a loudness rating of “most comfortable” is obtained. Complete real ear measurements at this final setting. Dorman et al. (2014)<sup>59</sup> found the loudness “balance point” by using a graphic of a head with an arc to indicate where sound is perceived. Adjustments were made until the sound was perceived at the centre.

**Note:** The appropriate adjustment/s—such as overall gain or compression—will be dependent on the patient’s

perception and the hearing aid settings that can be manipulated (device-specific). Best clinical judgment should be utilized.

**Caveat:** Caution should be taken if significant gain reduction is required as this may negatively impact audibility. The patient may need to be counseled regarding this and/or re-evaluation of cochlear implant settings may be recommended.



*The bimodal adjustments described in this paper can be carried out using the bimodal fitting panel in Genie. Use the loudness balance trimmer to increase or decrease the loudness across all frequencies simultaneously. Go to the Speech Rescue panel to apply frequency lowering. Use the High frequency On vs Off buttons to restrict the output bandwidth and/or give low frequency emphasis by sliding the low frequency emphasis trimmer.*

### Evaluation of benefit:

Documenting outcomes is important in order to demonstrate that at a minimum, performance in the bimodal condition is equivalent to or better than performance in the cochlear implant alone condition. There are many potential measures that may be employed.

**Recommended Strategy:** Word testing in quiet and sentence testing in quiet and in noise following the procedure outlined in the Minimum Speech Test Battery for Adults<sup>60</sup>. Age-appropriate measures should be utilized for younger children. If possible, evaluate in the CI-alone, HA-alone and Bimodal conditions. It should be noted that HA-alone performance may be very poor. On the contrary, if contralateral hearing could potentially contribute to speech perception performance, it is recommended that plugging and/or muffing that ear be considered to isolate the cochlear implant.

**Frequency Lowering:** If this is utilized, it may be helpful to also complete Ling six sound detection and discrimination/identification testing, particularly as it relates to /s/ and /S/.

**Loudness:** Consider aided soundfield thresholds for both the cochlear implant and the hearing aid. It should be emphasized that this is to be completed in addition to real ear verification. The purpose of utilizing this measure is to validate the detection levels of soft sounds and to compare minimum detection levels for soft inputs across the two devices. With a properly fitted cochlear implant, detection levels should be equivalent across the frequency range. With the hearing aid, detection levels will depend on the degree of hearing loss and hearing aid fitting. For example, high frequency thresholds are expected to be poor with a restricted bandwidth fitting. On the other hand, high frequency thresholds may be slightly better with frequency lowering versus wideband fitting (see Davidson et al., 2015). If the outcomes are not as expected, further evaluation of the fitting is recommended.

**Other Measures:** Current assessment measures focusing on speech recognition with standard clinical set-ups may not be sensitive enough to show significant bimodal benefit. If time permits and these measures are available, it may be worthwhile to consider the following:

- Speech perception testing with a multiple loudspeaker array (e.g. R-Space)
- Speech perception testing in reverberation<sup>61</sup>
- Sound localization
- Music Perception
- Speaker Recognition
- Emotion Recognition
- Subjective Questionnaires (e.g. Speech, Spatial, and Qualities of Hearing questionnaire, SSQ)

### What if no bimodal benefit is seen or bimodal decrement is observed?:

- Contact the CI audiologist to re-evaluate cochlear implant settings.
- Consider candidacy for bilateral cochlear implants.
- More extensive testing (e.g. for loudness balancing) may be necessary.
- More time may be needed to adapt to the bimodal signal<sup>62,63</sup> and/or bimodal fitting strategy (e.g. patient

used to wideband fitting is trying frequency lowering for the first time). Also, consider if aural rehabilitation would be beneficial if this is not already in place<sup>62</sup>.

• Reiss et al. (2014)<sup>64</sup> suggested that some bimodal users may experience significant difficulties with binaural spectral integration and are unable to resolve pitch mismatches perceived in the two ears. Alternative fitting strategies such as minimal-overlapping fitting<sup>65,66</sup> could be explored. Laria et al. (2014)<sup>67</sup> presented a case example wherein improvement was observed when the frequency range of the cochlear implant did not overlap with that of the hearing aid. It should be noted that current hybrid fitting protocols recommend some overlap between acoustic and electric hearing<sup>68,69</sup>. If this approach is appropriate for a particular bimodal patient, the audiologist may have to experiment to determine the edge frequency for acoustic and electric hearing that results in the best outcomes.

• Some recipients demonstrate bimodal decrement<sup>11,58</sup>. This is typically the exception, however, and not the rule.

### Concluding Remarks

Improvements in cochlear implant technology and outcomes have resulted in expanded candidacy criteria. As a result, more implant recipients present with significant residual hearing and bimodal stimulation should be considered in order to provide access to bilateral, binaural & complementary acoustic cues. At minimum, the hearing aid should be fit to target and verified using real ear measurements. Additional benefits may be obtained by further optimization of the hearing aid fitting. Optimization can include consideration of alternative frequency response settings and/or loudness balancing. Due to the inherent variability in the cochlear implant population, the approach that provides the most benefit will vary from patient to patient. Outcomes assessment is important in individualizing the hearing aid fitting approach.

### Case Studies:

Help us spread knowledge on bimodal fittings by adding your own bimodal fitting data to this paper. Send patient case description to Kamilla Angelo, PhD, [kian@oticon.com](mailto:kian@oticon.com).

- Patient 1 is an 88-year-old male with a long-standing progressive hearing loss.
- Patient 2 is a 57-year-old female with bilateral progressive hearing loss due to enlarged vestibular aqueducts.

## References

1. Dorman, M. F. & Gifford, R. H. Combining acoustic and electric stimulation in the service of speech recognition. *Int J Audiol* 49, 912-9 (2010).
2. Ching, T., Psarros, C., Hill, M., Dillon, H. & Incerti, P. Should Children Who Use Cochlear Implants Wear Hearing Aids in the Opposite Ear. *Ear & Hearing* 22, 365 -80 (2001).
3. Ching, T., Incerti, P & Hill, M. Binaural benefits for adults who use hearing aids and cochlear implants in opposite ears. *Ear and hearing* (2004). doi:10.1097/00003446-200402000-00002
4. Cullington, H. E. & Zeng, F.-G. G. Comparison of bimodal and bilateral cochlear implant users on speech recognition with competing talker, music perception, affective prosody discrimination, and talker identification. *Ear Hear* 32, 16-30 (2011).
5. Illg, A., Bojanowicz, M., Lesinski-Schiedat, A., Lenarz, T. & Büchner, A. Evaluation of the bimodal benefit in a large cohort of cochlear implant subjects using a contralateral hearing aid. *Otol. Neurotol.* 35, e240-4 (2014).
6. Chang, JE, Bai, JY & Zeng, FG. Unintelligible low-frequency sound enhances simulated cochlear-implant speech recognition in noise. *Biomedical Engineering* (2006). doi:10.1109/TBME.2006.883793
7. Kong, YY, Stickney, GS & Zeng, FG. Speech and melody recognition in binaurally combined acoustic and electric hearing. *The Journal of the Acoustical ...* (2005). doi:10.1121/1.1857526
8. Straatman, LV, Rietveld, A. & Beijen, J. Advantage of bimodal fitting in prosody perception for children using a cochlear implant and a hearing aid. *The Journal of the ...* (2010). at <<http://scitation.aip.org/content/asa/journal/jasa/128/4/10.1121/1.3474236>>
9. Sucher, CM & McDermott, HJ. Bimodal stimulation: benefits for music perception and sound quality. *Cochlear Implants International* (2009). doi:10.1002/cii.398
10. Hamzavi, J, Pok, M. S. & Gstoettner, W. Speech perception with a cochlear implant used in conjunction with a hearing aid in the opposite ear. .... *journal of audiology* (2004). doi:10.1080/14992020400050010
11. Tyler, RS, Parkinson, AJ, Wilson, BS & Witt, S. Patients utilizing a hearing aid and a cochlear implant: speech perception and localization. *Ear and ...* (2002). doi:10.1097/00003446-200204000-00003
12. Armstrong, M., Pegg, P., James, C. & Blamey, P. Speech perception in noise with implant and hearing aid. *The American journal of otology* 18, S140-1 (1997).
13. Blamey, P.J., Armstrong, M. & James, J. Cochlear implants, hearing aids, or both together? In G. M. Clark (Ed.). *Cochlear implants* 273 - 277 (1997).
14. Berrettini, S., Passetti, S., Giannarelli, M. & Forli, F. Benefit from bimodal hearing in a group of prelingually deafened adult cochlear implant users. *American journal of otolaryngology* 31, 332-338 (2010).
15. Ullauri, A, Crofts, H & Wilson, K. Bimodal benefits of cochlear implant and hearing aid (on the non-implanted ear): a pilot study to develop a protocol and a test battery. *Cochlear implants ...* (2007). doi:10.1179/cim.2007.8.1.29
16. Nittrouer, S, Caldwell, A, Lowenstein, JH & Tarr, E. Emergent literacy in kindergartners with cochlear implants. *Ear and ...* (2012). at <<http://www.ncbi.nlm.nih.gov/pmc/articles/PMC3419773/>>
17. Nittrouer, S & Chapman, C. The effects of bilateral electric and bimodal electric-Acoustic stimulation on language development. *Trends in amplification* (2009). at <<http://tia.sagepub.com/content/13/3/190.short>>
18. Nittrouer, S., Kuess, J. & Lowenstein, J. H. Speech perception of sine-wave signals by children with cochlear implants. *J. Acoust. Soc. Am.* 137, 2811 (2015).
19. Offeciers, E, Morera, C, Müller, J & Huarte, A. International consensus on bilateral cochlear implants and bimodal stimulation. *Acta oto- ...* (2005). at <<http://www.ncbi.nlm.nih.gov/pubmed/16109670>>
20. Balkany, T. et al. William House Cochlear Implant Study Group: position statement on bilateral cochlear implantation. *Otology & Neurotology* 29, 107-108 (2008).
21. Scherf, F. & Arnold, LP. Exploring the clinical approach to the bimodal fitting of hearing aids and cochlear implants: results of an international survey. *Acta oto-laryngologica* (2014). at <<http://informahealthcare.com/doi/abs/10.3109/00016489.2014.914244>>
22. Siburt, HW & Holmes, AE. Bimodal Fitting: A survey of current clinical practice. *American journal of audiology* (2015). at <<http://aja.pubs.asha.org/article.aspx?articleid=2279994>>
23. Keilmann, AM, Bohnert, AM & Gosepath, J. Cochlear implant and hearing aid: a new approach to optimizing the fitting in this bimodal situation. *European Archives of Oto- ...* (2009). doi:10.1007/s00405-009-0993-9
24. Huart, SA & Sammeth, CA. Hearing aids plus cochlear implants: Optimizing the bimodal paediatric fitting. *The Hearing Journal* (2008). doi:10.1097/01.HJ.0000342441.45181.6e
25. Yehudai, N, Shpak, T, Most, T & Luntz, M. Functional status of hearing aids in bilateral-bimodal users. *Otology & Neurotology* (2013). doi:10.1097/MAO.0b013e3182898131
26. Harris, MS & Hay McCutcheon, M. An analysis of hearing aid fittings in adults using cochlear implants and contralateral hearing aids. *The Laryngoscope* (2010). doi:10.1002/lary.21148
27. American Academy of Audiology. Guidelines for the Audiologic Management of Hearing Impairment. . *Audiology Today* 18, (2006).
28. American Academy of Audiology. Clinical Practice Guidelines: Paediatric Amplification. . [http://audiology-web.s3.amazonaws.com/migrated/PaediatricAmplificationGuidelines.pdf\\_539975b3e7e9f1.74471798.pdf](http://audiology-web.s3.amazonaws.com/migrated/PaediatricAmplificationGuidelines.pdf_539975b3e7e9f1.74471798.pdf)
29. Ching, T. Y., van Wanrooy, E. & Dillon, H. Binaural-bimodal fitting or bilateral implantation for managing severe to profound deafness: a review. *Trends Amplif* 11, 161-92 (2007).
30. Ching, T. Y., Incerti, P. & Plant, K. Electric-acoustic stimulation: for whom, in which ear, and how. *Cochlear Implants Int* 16 Suppl 1, S12-5 (2015).
31. Morera, C, Cavalle, L & Manrique, M. Contralateral hearing aid use in cochlear implanted patients: multicenter study of bimodal benefit. *Acta oto- ...* (2012). at <<http://informahealthcare.com/doi/abs/10.3109/00016489.2012.677546>>
32. Davidson, LS, Firszt, JB & Brenner, C. Evaluation of Hearing Aid Frequency Response Fittings in Paediatric and Young Adult Bimodal Recipients. *Journal of the ...* (2015). at <<http://www.ingentaconnect.com/content/aaa/jaaa/2015/00000026/00000004/art00007>>
33. Stelmachowicz, PG & Pittman, AL. The importance of high-frequency audibility in the speech and language development of children with hearing loss. ... -*Head & Neck ...* (2004). doi:10.1001/archotol.130.5.556
34. Seeber, BU, Baumann, U & Fastl, H. Localization ability with bimodal hearing aids and bilateral cochlear implants. *The Journal of the Acoustical ...* (2004). doi:10.1121/1.1776192
35. Potts, LG, Skinner, MW & Litovsky, RA. Recognition and localization of speech by adult cochlear implant recipients wearing a digital hearing aid in the nonimplanted ear (bimodal hearing). *Journal of the ...* (2009). at <<http://www.ncbi.nlm.nih.gov/pmc/articles/PMC2876351/>>
36. Bartov, T & Most, T. Song Recognition by Young Children With Cochlear Implants: Comparison Between Unilateral, Bilateral, and Bimodal Users. *Journal of Speech* (2014). at <<http://sig17perspectives.pubs.asha.org/article.aspx?articleid=1870320>>
37. Fata, E. F., James, CJ & Laborde, ML. How much residual hearing is 'useful' for music perception with cochlear implants? *Audiology and ...* (2009). doi:10.1159/000206491
38. Gfeller, K, Oleson, J, Knutson, JF & Breheny, P. Multivariate predictors of music perception and appraisal by adult cochlear implant users. *Journal of the ...* (2008). at <<http://www.ncbi.nlm.nih.gov/pmc/articles/PMC2677551/>>
39. Neuman, AC & Svirsky, MA. Effect of hearing aid bandwidth on speech recognition performance of listeners using a cochlear implant and contralateral hearing aid (bimodal hearing). *Ear and hearing* (2013). at <<http://europepmc.org/articles/pmc3748228>>
40. Sheffield, SW & Gifford, RH. The benefits of bimodal hearing: effect of frequency region and acoustic bandwidth. *Audiology and Neurotology* (2014). doi:10.1159/000357588


41. Kong, YY & Carlyon, RP. Improved speech recognition in noise in simulated binaurally combined acoustic and electric stimulations). *The Journal of the Acoustical Society of ...* (2007). doi:10.1121/1.2717408
42. Zhang, T, Dorman, MF & Spahr, AJ. Information from the voice fundamental frequency (F0) region accounts for the majority of the benefit when acoustic stimulation is added to electric stimulation. *Ear and hearing* (2010). at <<http://www.ncbi.nlm.nih.gov/pmc/articles/PMC3684557/>>
43. Shpak, T, Most, T & Luntz, M. Fundamental frequency information for speech recognition via bimodal stimulation: cochlear implant in one ear and hearing aid in the other. *Ear and hearing* (2014). doi:10.1097/AUD.0b013e3182a2c814
44. Zhang, T., Dorman, M. F., Gifford, R. & Moore, B. C. Cochlear dead regions constrain the benefit of combining acoustic stimulation with electric stimulation. *Ear Hear* 35, 410-7 (2014).
45. Mok, M., Grayden, D., Dowell, R. C. & Lawrence, D. Speech perception for adults who use hearing aids in conjunction with cochlear implants in opposite ears. *J. Speech Lang. Hear. Res.* 49, 338-51 (2006).
46. Mok, M, Galvin, KL, Dowell, RC & McKay, CM. Speech perception benefit for children with a cochlear implant and a hearing aid in opposite ears and children with bilateral cochlear implants. *Audiology and Neurotology* (2010). doi:10.1159/000219487
47. Preminger, JE & Carpenter, R. A clinical perspective on cochlear dead regions: intelligibility of speech and subjective hearing aid benefit. *Journal of the American ...* (2005). at <<http://www.ingentaconnect.com/content/aaa/jaaa/2005/00000016/00000008/art00009>>
48. Glista, D, Scollie, S, Bagatto, M & Seewald, R. Evaluation of nonlinear frequency compression: Clinical outcomes. .... *Journal of Audiology* (2009). doi:10.1080/14992020902971349
49. McCreery, RW, Brennan, MA, Hoover, B & Kopun, J. Maximizing audibility and speech recognition with non-linear frequency compression by estimating audible bandwidth. *Ear and ...* (2013). at <<http://www.ncbi.nlm.nih.gov/pmc/articles/PMC3566286/>>
50. Hua, H, Johansson, B & Jönsson, R. Cochlear implant combined with a linear frequency transposing hearing aid. *Journal of the ...* (2012). at <<http://www.ingentaconnect.com/content/aaa/jaaa/2012/00000023/00000009/art00006>>
51. Perreau, A. E., Bentler, R. A. & Tyler, R. S. The contribution of a frequency-compression hearing aid to contralateral cochlear implant performance. *J Am Acad Audiol* 24, 105-20 (2013).
52. McDermott, H & Henshall, K. The use of frequency compression by cochlear implant recipients with postoperative acoustic hearing. *Journal of the American Academy ...* (2010). at <<http://www.ingentaconnect.com/content/aaa/jaaa/2010/00000021/00000006/art00005>>
53. Park, LR, Teagle, H., Buss, E & Roush, PA. Effects of frequency compression hearing aids for unilaterally implanted children with acoustically amplified residual hearing in the nonimplanted ear. *Ear and ...* (2012). doi:10.1097/AUD.0b013e31824a3b97
54. Alexander, J. M. 20Q: The Highs and lows of frequency lowering amplification. . *AudiologyOnline Article* 11772, (2013).
55. Blamey, P. J., Dooley, G. J., James, C. J. & Parisi, E. S. Monaural and binaural loudness measures in cochlear implant users with contralateral residual hearing. *Ear Hear* 21, 6-17 (2000).
56. Van Hoesel, R. Contrasting benefits from contralateral implants and hearing aids in cochlear implant users. *Hearing research* (2012). at <<http://www.sciencedirect.com/science/article/pii/S0378595511002942>>

57. Francart, T & McDermott, HJ. Psychophysics, fitting, and signal processing for combined hearing aid and cochlear implant stimulation. *Ear and hearing* (2013). doi:10.1097/AUD.0b013e31829d14cb
58. Litovsky, RY & Johnstone, PM. Benefits of bilateral cochlear implants and/or hearing aids in children: Beneficios de los implantes cocleares bilaterales y/o auxiliares auditivos en niños. *International Journal of ...* (2006). doi:10.1080/14992020600782956
59. Dorman, MF, Loizou, P, Wang, S & Zhang, T. Bimodal cochlear implants: the role of acoustic signal level in determining speech perception benefit. *Audiology and ...* (2014). doi:10.1159/000360070
60. Fabry, D., Firszt, J. B., Gifford, R. H., Holden, L. K. & Koch, D. B. Evaluating speech perception benefit in adult cochlear implant recipients. *Audiology Today May*, 36 - 43 (2009).
61. Mason, M & Kokkinakis, K. Perception of consonants in reverberation and noise by adults fitted with bimodal devices. *Journal of Speech* (2014). at <<http://jshd.pubs.asha.org/article.aspx?articleid=1832502>>
62. Zhang, T., Dorman, M. F., Fu, Q.-J. & Spahr, A. J. Auditory training in patients with unilateral cochlear implant and contralateral acoustic stimulation. *Ear Hear* 33, e70-9 (2012).
63. Luntz, M, Shpak, T & Weiss, H. Binaural-bimodal hearing: concomitant use of a unilateral cochlear implant and a contralateral hearing aid. *Acta oto-laryngologica* (2005). at <<http://informahealthcare.com/doi/abs/10.1080/00016480510035395>>
64. Reiss, L., Ito, RA, Eggleston, JL & Wozny, DR. Abnormal binaural spectral integration in cochlear implant users. *Journal of the Association for ...* (2014). doi:10.1007/s10162-013-0434-8
65. Green, T, Faulkner, A & Rosen, S. Overlapping frequency coverage and simulated spatial cue effects on bimodal (electrical and acoustical) sentence recognition in noise. .... *Journal of the Acoustical Society of ...* (2014). doi:10.1121/1.4861843
66. Reiss, L., Ito, RA, Eggleston, JL & Liao, S. Pitch Adaptation Patterns in Bimodal Cochlear Implant Users: Over Time and After Experience. *Ear and ...* (2015). doi:10.1097/AUD.000000000000114
67. Laria, C, Auletta, G, Riccardi, P & Papa, C. ... good performance with bimodal stimulation in a like-hybrid modality in a patient with profound bilateral sensorineural hearing loss with low-frequencies preservation. *American journal of ...* (2014). at <<http://www.sciencedirect.com/science/article/pii/S019607091300207X>>
68. Gifford, R. 20Q: Hybrid/EAS cochlear implants- New research and clinical tips. *AudiologyOnline Article* 828, (2011).
69. Vermeire, K, Anderson, I & Flynn, M. The influence of different speech processor and hearing aid settings on speech perception outcomes in electric acoustic stimulation patients. *Ear and ...* (2008). doi:10.1097/AUD.0b013e31815d6326

**Case Studies: Patient 1**

Adult Case Study: Patient 1 is an 88-year-old male with a long-standing progressive hearing loss. He has a history of occupational and military noise exposure without the use of hearing protection. At the time he went through the cochlear implant candidacy process, he had normal hearing at 125-250 Hz steeply sloping to a profound sensorineural hearing loss in the higher frequencies (See Figure 1). He had been utilizing hearing aids for over 20 years. He met the candidacy criteria for cochlear implantation and received a cochlear implant on his left ear.

Figure 1. Pre-implant audiogram for Patient 1.

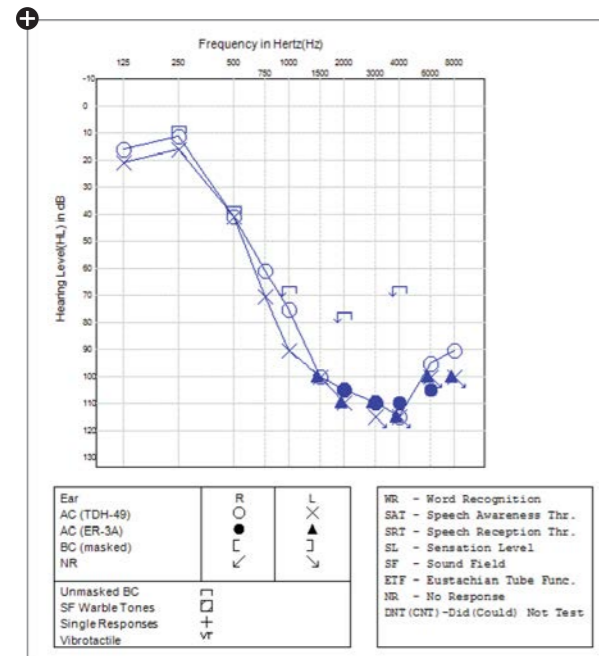

Patient 1 continues to utilize a hearing aid in the non-implanted ear. This currently utilized hearing aid has frequency lowering (transposition) enabled to maximize high frequency audibility. The patient has indicated a preference for frequency lowering versus standard fitting. Probe microphone measurements indicate audibility up to about 1000 Hz (Figure 2). Using the Frequency Lowering Fitting Assistant (Figure 3), best performance was predicted using a start frequency (the first frequency that will be transposed) of 1600 Hz and “Expanded Mode” (transposition of 5 frequency bands versus 3 bands in “Basic Mode”). Gain for transposed sound is set at a default of 0.

Figure 2. Speech Map showing the probe microphone real ear measurement for Patient 1’s Hearing Aid. Test 1 (green) shows results for standard wideband fitting while Test 4 (orange) shows results for frequency lowering.


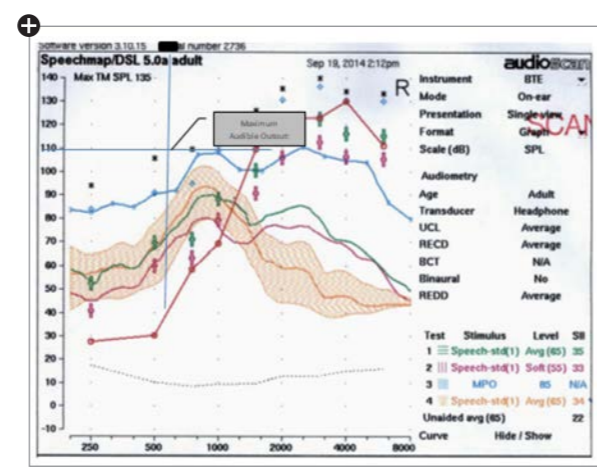


Figure 3. Settings provided by the Frequency Lowering Fitting Assistant.

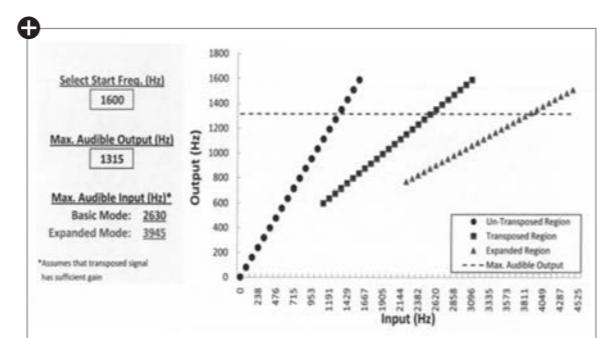


Patient 1 has had his implant for about a year and speech perception performance is shown in Table 1. Results indicate improvement over pre-implant performance and best performance in the bimodal condition. The patient is pleased with his current hearing.

TEST	Pre-Implant Score for Left Hearing Aid	Pre-Implant Score for Both Hearing Aids	Score with the Right Hearing Aid Alone	Score with Cochlear Implant Alone	Score in the Bimodal Condition
CNC Test					
-Words Correct	12%	Did not Test	24%	40%	40%
-Phonemes Correct	31%		44%	53%	71%
HINT Sentences					
-Words Correct in Quiet	26%	62%	43%	57%	92%
AZ Bio Sentences					
-Words Correct in Quiet	Did not Test	32%	Did not Test	27%	67%


**Case Studies: Patient 2**

Patient 2 is a 57-year-old female with bilateral progressive hearing loss due to enlarged vestibular aqueducts. At the time she went through the cochlear implant candidacy process she had a mild sloping to profound hearing loss in the right ear and a severe to profound loss in the left ear. She met the candidacy criteria for cochlear implantation and received an implant in her left ear. She has had her device for approximately a year and a half.

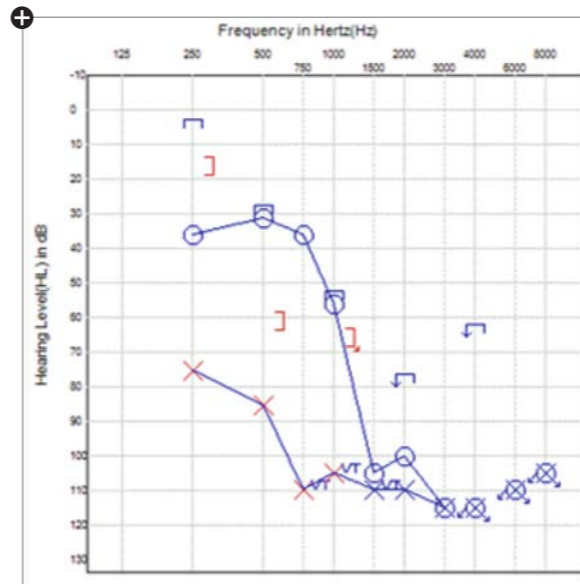

Patient 2 continued utilizing her hearing aid following cochlear implantation. Her hearing aid is fitted to match adult DSL5.0 targets for as wide a bandwidth as possible given her steeply sloping loss. No frequency lowering is enabled even though it is available in her instrument because NFC will not improve audibility since MAF is 1500 Hz and 1500 Hz is the lowest start frequency available.

Due to binaural loudness summation, the patient requested a slight gain decrease. Probe microphone real ear measurement results are shown below.

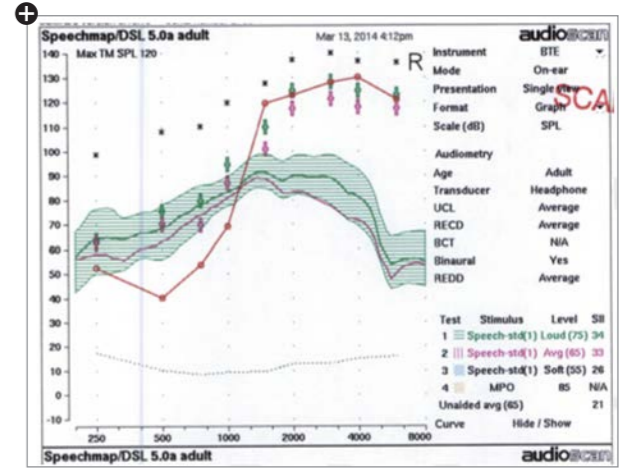

Speech perception test results are shown in Table 2. Results indicate significant improvement in scores over pre-implant performance. Performance in the bimodal condition is similar to or slightly improved over performance with the cochlear implant alone. The patient is satisfied with her current hearing.

TEST	Pre-Implant Score for Left Hearing Aid	Pre-Implant Score for Both Hearing Aids	Score with the Right Hearing Aid Alone	Score with Cochlear Implant Alone	Score in the Bimodal Condition
CNC Test					
-Words Correct	Did not Test	38%	Did not Test	80%	85%
-Phonemes Correct		80%		91%	95%
HINT Sentences					
-Words Correct in Quiet	32%	51%	58%	90%	98%
AZ Bio Sentences					
-Words Correct in Quiet	Did not Test	26%	Did not Test	90%	98%
-Words Correct in Noise		Did not Test		66%	76%



[oticon.global](http://oticon.global)

**oticon**  
PEOPLE FIRST